


# Defeating Pain with MR4 Laser

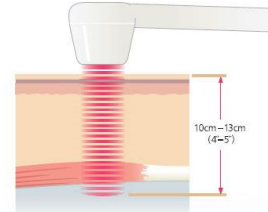

**Pulse drives the light energy, to the target tissue, up to 10-13cm (4-5 inches) deep**

## Revolutionary MR4 Super Pulsed Laser Technology for treatment of pain and injuries

### Advanced Non-invasive Pain Relief Device

- Anti-inflammatory action
- Increases vascular activity/metabolic activity
- Stimulates nerve function
- Rapid cell growth
- Reduces fibrous tissue formation

### Indication includes:

- Acute, sub-acute and chronic pain
- Shoulder pain
- Carpal tunnel
- Back pain, Tennis elbow
- Arthritic pain
- Many other indications

### 4 Radiances Technology

Combined into a **SINGLE** unit

- **The GaAs Super Pulse Laser / 905nm**
- **Four infrared diodes / 875nm**
- **Four red diodes / 660nm**
- **Static Magnetic Field/enhance energy potential at the molecular levels**

### **Super Pulsed Laser Power**

Radiation Penetrates  
Peak Pulse Power 25,000mW  
or 50,000mW in Laser Shower (20 cm<sup>2</sup>)  
10-13cm (4-5 inches) deep

**Also Great for Anti-Aging and Cosmetic Solutions**

Item #	Product Description	Professional
17-9211	MR4 Laser (SE25 – 25W)	<b>\$8,695</b>
17-9212	MR4 LS50 Shower Laser (4 x 905nm diodes)	<b>\$9,950</b>
17-9213	MR4 Laser Stim with Target Technology	<b>\$11,500</b>
17-9214	MR4 ULTRA Package Laser (Stim+Shower)	<b>\$14,500</b>
17-9210	TerraQuant TQ Solo (15W - Portable)	<b>\$3,750</b>
17-9215	MR4 ACTIV Portable Cordless Laser (25W)	<b>\$9,950</b>

### Advantages of MR4 and LaserStim™ with TARGET™

- Synergistic effect laser therapy provides benefits of light therapy and neurostimulation, identifying optimal treatment sites, and provides highly targeted treatment where photons need to be delivered and electro-stimulation needs to be applied.
- Faster pain relief and reduction in recovery time
- Constantly changing impulses, preventing accommodation of tissue to neurostimulation
- Deepest tissue penetration of photonic energy due to Super Pulsed Laser
- Preset frequencies allow clinician to control depth of therapeutic energy from the surface up to 13 cm (5 inches)
- Can also target acupuncture points

### Technological Advantage: LaserSweep™

The continual use of the same frequency, as found in continuous wave lasers, can lead ultimately to the beginning of adaptation processes and thereby can diminish the therapeutic effect. The Multi Radiance MR4™ is capable of creating unique and exclusive, patented "sweep" laser programs.

With each laser impulse of a Sweep Program, the "frequency" is changed, varying between end points of pre-selected frequencies. This effectively change the laser's target depth of penetration, in essence, "scanning" the target tissue from superficial to deep, while preventing biological adaptation to the treatment.

# Kosim

ideal for health

**KOSIM Group Company Inc.**

55 Valleywood Drive, Markham, Ontario L3R 5L9 Canada

T: 905-475-2380 F: 905-475-2303

[www.kosim.ca](http://www.kosim.ca)

Email: [sales@kosim.ca](mailto:sales@kosim.ca)

Toll-Free: 1-877-50-KOSIM

**WE SUPPLY A RANGE OF CLINICAL SUPPLIES & EQUIPMENT**

## ACCESSORIES FOR YOUR MR4



### 17-9210 Portable TQ Solo

**\$3,750**

Laser Peak Power: 15W (15,000mw)

Infrared LED's: 60mW

Red LED's: 7.5mW

Static Magnetic Field

**Wavelength:** Laser 905nm,

Infrared 875nm, & Red 660nm

Class II Equipment Applied Part Type BF

Laser Class: 1M (IEC 60825-1)

Weight: 300g , Dimensions: 5.5 x 2 x 2"



Most versatile, perfect for use with probes

Provides complete spectral irradiance

### 17-9202 Extra SE25 Standard (TQ1) Emitter

**\$2,895**

Laser Peak Power: 25W (25,000mw)

Infrared LED's: 60mW

Red LED's: 7.5mW

Static Magnetic Field

Coverage: 4 cm<sup>2</sup>

**Wavelength:**

1 x Laser 905nm, 4 x Infrared 875nm, & 4 x Red 660nm

Class II Equipment Applied Part Type BF



Covers 5 times more area to treat large muscle groups

Incorporates 4 super pulsed GaAs diodes

### 17-9203 Extra LS50 Laser Shower Emitter

**\$4,150**

Laser Peak Power: 50W (50,000mw)

Infrared LED's: 60mW

Red LED's: 7.5mW

Static Magnetic Field

Coverage: 20 cm<sup>2</sup>

**Wavelength:**

4 x Laser 905nm, 4 x Infrared 875nm, & 4 x Red 660nm

Class II Equipment Applied Part Type BF



Highly effective and reimbursable

Combines light therapy and electrical stimulation

### 17-9209 Extra LaserStim Emitter

**\$5,495**

Laser Peak Power: 25W (25,000mw)

Infrared LED's: 60mW

Red LED's: 7.5mW

Static Magnetic Field + Electrical Stimulation

Coverage: 4 cm<sup>2</sup> (ATTENDED THERAPY ONLY)

**Wavelength:**

1 x Laser 905nm, 4 x Infrared 875nm, & 4 x Red 660nm

Class II Equipment Applied Part Type BF



### 17-9201 Set of 4 Acupuncture Point Probe

**\$540**

The probes allow the laser cluster to be adapted to treat smaller target areas without having to purchase additional expensive equipment.



### 17-9208 MR4 Deluxe Magna Cart with Double Arms **\$1,680**

Your choice of hands-free armatures allows for unattended treatments and is easily attachable to any cart or table. Each armature is designed to provide the doctor with a hands-free environment while working with Multi Radiance Medical Laser emitters.

Instead of doctors holding the emitter on the treatment area, the hands-free armature will do it for them. An optional C-clamp allows easy attachment to Multi Radiance Medical carts and trolleys.

- Size: 34" Height
- Accessories tray include
- Lightweight and easy to use
- Compatible with SE25 and LS50 emitters